



P R O S T H E T I C S

SEATTLE SELECT
STANCE FLEXION
KNEE

4-HOLE
ROTATABLE
FLANGE



AAA
DOUBLE-ENDED
ADAPTERS



ALUMINUM PYLON
BONDED to
MALE ADAPTER

SEATTLE
TRIUMPH



Delaney
Varsity Cheerleader,
Soccer Team member

From a miracle baby to an active, energetic athlete and young woman! Delaney, born premature and without her left foot, could never be stopped. She grew up wearing Seattle Child's Play Feet, and now as a young woman, she has grown into wearing a Trulife adult Seattle Foot. Delaney is now attending a university to study early childhood development and will also be teaching at a cheerleading camp.

Contents

FEET

Seattle Triumph	4
Seattle LP	5
Seattle Zumo 450	6
Seattle Energy	7
Seattle Lightfoot with Pyramid	8
Seattle Kinetic Edge with Pylon	9
Seattle Kinetic Light with Pylon	10
Seattle Kinetic with Pylon	11
Seattle Natural with Pyramid	12
Seattle SACH with Pyramid	13

KNEES

Seattle Select Stance Flexion Knee	14
Seattle Natural Stance Flexion Knee	15
Seattle PC Friction Control Knee with Lock	15

AAA MODULAR COMPONENTS

3-Prong Adapters	16
4-Hole Adapters	16
Clamp Adapters	16
Double-Ended Adapters	17
Foot Pyramid	17
Seattle Mounting Plate	17
Pylons	17
Heavy-Duty Titanium Adapters	18

OTHER

Green Dot Valves	19
Below Knee Valve	19
Seattle Endoskeletal Cosmesis	20
Seattle Ankle Block	20
Seattle Symes Nut	21
Seattle Mounting Plate	21

CHILD'S PLAY

Seattle Child's Play Energy	22
Seattle Child's Play Knee with Manual Lock	23
Child's Play Adapters	24
Child's Play Endo/Exo Components	25

WARRANTY

Warranty	27
----------	----



Seattle Triumph

- Narrowed shank provides multi-axial compliance to uneven terrain increasing stability and limb comfort
- 9 categories to choose from
- Optional toe wedges provide increased dorsiflexion resistance and improved proprioception at mid-stance for increased stability
- Multi-axial rotation and shock absorption
- Full-length split keel
- Smooth rollover and controlled plantar flexion
- Pre-assembled with Spectra sock in removable foot shell
- Foot shell available in light and dark

Indications: Medium to high impact levels, K3-K4, maximum weight of 160 kg /350 lb*

Selection

Model No.	Description	Size	Weight Limit
STF190	Seattle Triumph, Light Foot Shell	22-30 cm	160 kg / 350 lb*
STF193	Seattle Triumph, Dark Foot Shell	22-30 cm	160 kg / 350 lb*
SFC290	Replacement Foot Shell, Light	22-30 cm	N/A
SFC293	Replacement Foot Shell, Dark	22-30 cm	N/A

Specifications

Weight Limit	160 kg / 350 lb*
Product Sizes	22-30 cm
Toe	Solid Toe
Foot Shell	Light, Dark
Product Weight	493 g / 1 lb (Size 25 without Foot Shell)
Min. Clearance	18 cm / 7.1 in
Heel Rise	9.5 mm / 3/8 in
Warranty Keel: Foot Shell:	3 years 6 months
Pre-assembled with	Foot Shell and Spectra Sock
Made in USA	

*For low impact level, category 9

Category Selection Instructions

To optimize the selection and ensure amputee's safety, follow the two steps below to determine the appropriate category.

1. Locate the column that corresponds with the amputee's impact level.
2. Within the selected column, locate the amputee's weight.

Important Note: If the amputee has a long BK, carries heavy loads or will reach a higher impact level within a year, choose the next category higher.

Warning: Choosing a lower strength category than what is suggested based on the above procedure and patient data will void the warranty. If your patient's weight exceeds the limits of the chart please call Trulife Customer Service.

Category	9	8	7	6	5	4	3	2	1
Low Impact Walking, Uneven Surfaces	141-160 kg 310-350 lb	121-140 kg 266-309 lb	101-120 kg 221-265 lb	91-100 kg 199-220 lb	81-90 kg 177-198 lb	71-80 kg 155-176 lb	66-70 kg 144-154 lb	60-65 kg 132-143 lb	< 60 kg < 132 lb
Medium Impact Light Sports	124-145 kg 271-320 lb	106-123 kg 232-270 lb	89-105 kg 194-231 lb	80-88 kg 175-193 lb	71-79 kg 155-174 lb	62-70 kg 136-154 lb	58-61 kg 126-135 lb	53-57 kg 117-125 lb	< 53 kg < 117 lb
High Impact Running, Basketball	106-125 kg 232-274 lb	91-105 kg 199-231 lb	76-90 kg 166-198 lb	69-75 kg 150-165 lb	61-68 kg 133-149 lb	53-60 kg 117-132 lb	50-53 kg 108-116 lb	46-49 kg 101-107 lb	< 46 kg < 101 lb



AAA100
Adjustable Height
Clamp Adapter
(page 16)

AAA135
Variable Height
Adapter
(page 17)



Seattle LP

- Low profile foot for long residual limb
- Multi-axial
- C-shape replicates true, anatomical plantar flexion
- Full length, split, composite keel for stability on uneven surfaces
- Removable foot shell

Indications: Medium to high impact level, K3-K4, maximum weight of 136 kg / 300 lb*; ideal for amputees with long residual limbs

Selection

Model No.	Description	Size	Weight Limit
SLP350	Seattle LP, Light Foot Shell	22-30 cm	136 kg / 300 lb*
SLP353	Seattle LP, Dark Foot Shell	22-30 cm	136 kg / 300 lb*
SFC290	Replacement Foot Shell, Light	22-30 cm	N/A
SFC293	Replacement Foot Shell, Dark	22-30 cm	N/A

Specifications

Weight Limit	136 kg / 300 lb*
Product Sizes	
Keel 1-3:	22-30 cm
Keel 4:	26-30 cm
Toe	Solid Toe
Foot Shell	Light, Dark
Product Weight	400 g / 14 oz (Size 25 without Foot Shell)
Min. Clearance	9 cm / 3.6 in
Heel Rise	9.5 cm / 3/8 in
Warranty	
Keel:	3 years
Foot Shell:	6 months
Includes	Foot Shell, Spectra Sock, Bearing Kit
Made in USA	

*For medium impact level, category 4

Category Selection Instructions

To optimize the selection and ensure amputee's safety, follow the two steps below to determine the appropriate category.

1. Locate the column that corresponds with the amputee's impact level.
2. Within the selected column, locate the amputee's weight.

Important Note: If the amputee has a long BK, carries heavy loads or will reach a higher impact level within a year, choose the next category higher.

Warning: Choosing a lower strength category than what is suggested based on the above procedure and patient data will void the warranty. If your patient's weight exceeds the limits of the chart, please browse the catalog for other foot options or call Trulife Customer Service.

	Low-Medium	Medium-High	High
Category	Walking Uneven Surfaces	Light Sports	Running, Basketball
4 (Sizes 26-30)	101-136 kg 221-300 lb	88-118 kg 193-260 lb	81-109 kg 179-240 lb
3 (Sizes 22-30)	81-100 kg 177-220 lb	74-87 kg 162-192 lb	68-80 kg 150-176 lb
2 (Sizes 22-30)	61-80 kg 133-176 lb	56-73 kg 122-161 lb	51-67 kg 111-149 lb
1 (Sizes 22-30)	< 60 kg < 132 lb	< 55 kg < 121 lb	< 50 kg < 110 lb



AAA100
Adjustable Height
Clamp Adapter
(page 16)

AAA135
Variable Height
Adapter
(page 17)



Seattle Zumo 450

- Heavy-duty foot withstands up to 450 lb
- Durable composite design for K3 users
- Pre-assembled with a male pyramid for quick and easy installation
- Split composite keel for increased stability
- Dynamic toe and cushioned heel resposns
- Removable, life-like sandal-toe foot shell

Indications: Medium impact level, K3, maximum weight of 204 kg / 450 lb

Selection

Model No.	Description	Size	Weight Limit
SZM450	 Seattle Zumo, Light Foot Shell	22-30 cm	204 kg / 450 lb
SZM453	 Seattle Zumo, Dark Foot Shell	22-30 cm	204 kg / 450 lb
SFC190	 Replacement Foot Shell, Light	22-30 cm	N/A
SFC193	 Replacement Foot Shell, Dark	22-30 cm	N/A

Specifications

Weight Limit	204 kg / 450 lb
Product Sizes	22-30 cm
Toe	Sandal-Toe
Foot Shell	Light, Dark
Product Weight	750 g / 26.5 oz (Size 27 with Foot Shell)
Min. Clearance	
Heel Rise	9.5 cm / 3/8 in
Warranty Keel:	2 years
Foot Shell:	6 months
Includes	Foot Shell, Spectra Sock
Made in USA	



Trulife Complementary Items



AAAHTi225
Heavy-Duty Adjustable
Height Clamp Adapter
(page 18)



AAAHTi200
Heavy-Duty Bonded
Pylon, 16 inch
(page 18)



Seattle Energy

- Dynamic toe and cushioned heel response
- Full length, split composite keel for increased stability
- Removable, life-like, sandal-toe foot shell
- Pre-assembled with male pyramid for easy installation
- Pyramid and spacer can be removed for Symes patient application

Indications: Medium impact level, K3, maximum weight of 136 kg / 300 lb*

Selection

Model No.	Description	Size	Weight Limit
SEF190	Seattle Energy, Light Foot Shell	22-30 cm	136 kg / 300 lb*
SEF193	Seattle Energy, Dark Foot Shell	22-30 cm	136 kg / 300 lb*
SFC190	Replacement Foot Shell, Light	22-30 cm	N/A
SFC193	Replacement Foot Shell, Dark	22-30 cm	N/A

Specifications

Weight Limit	136 kg / 300 lb*
Product Sizes	22-30 cm
Toe	Sandal Toe
Foot Shell	Light, Dark
Product Weight	330 g / 0.7 lb (Size 25 without Foot Shell)
Min. Clearance	10 cm / 4 in
Heel Rise	9.5 mm / 3/8 in
Warranty Keel:	2 years
Foot Shell:	6 months
Pre-assembled with	Foot Shell, Spectra Sock, Foot Bolt, 2.5 cm Foot Spacer and Foot Pyramid
Made in USA	

*For medium impact level, category 4

Category Selection Instructions

To optimize the selection and ensure amputee's safety, select the appropriate category based on the amputee's weight. Chart accommodates low-medium impact level amputees.

Important Note: If the amputee has a long BK, carries heavy loads or will reach a higher impact level within a year, choose the next category higher.

Warning: Choosing a lower strength category than what is suggested based on the above procedure and patient data will void the warranty. If your patient's weight exceeds the limits of the chart please browse the catalog for other foot options or call Trulife Customer Service.

	Low	Medium-Low	Medium
Category	Walking with Aid	Limited Walking	Walking Uneven Surfaces
4	101-136 kg 221-300 lb	101-136 kg 221-300 lb	101-136 kg 221-300 lb
3	81-100 kg 177-220 lb	81-100 kg 177-220 lb	81-100 kg 177-220 lb
2	61-80 kg 133-176 lb	61-80 kg 133-176 lb	61-80 kg 133-176 lb
1	< 60 kg < 132 lb	< 60 kg < 132 lb	< 60 kg < 132 lb

Trulife Complementary Items



AAA100
Adjustable Height
Clamp Adapter
(page 16)

AAA135
Variable Height
Adapter
(page 17)

Seattle Lightfoot² with Pyramid

- Arrives ready to install
- Natural looking, sandal toe foot shell (non-removable)
- Pyramid and spacer can be removed for Symes patient application

Indications: Medium impact level, K3, maximum weight of 136 kg / 300 lb*

Selection

Model No.	Description	Size	Weight Limit
SLF195	Seattle Lightfoot ² with Pyramid 9.5 mm / 3/8 in Heel, Light Foot Shell	22-30 cm	136 kg / 300 lb*
SLF198	Seattle Lightfoot ² with Pyramid 9.5 mm / 3/8 in Heel, Dark Foot Shell	22-30 cm	136 kg / 300 lb*



Specifications

Weight Limit	136 kg / 300 lb*
Product Sizes	22-30 cm
Toe	Sandal Toe
Foot Shell	Light, Dark
Product Weight	480 g / 1 lb (Size 25)
Min. Clearance	9 cm / 3.5 in
Heel Rise	9.5 mm / 3/8 in
Warranty	
Keel:	2 years
Foot Shell:	1 year
Pre-assembled with	Foot Bolt, 2.5 cm Foot Spacer Foot Pyramid
Made in USA	

*For medium impact level

Category Selection Instructions

To optimize the selection and ensure amputee's safety, follow the two steps below to determine the appropriate keel.

1. Locate the column that corresponds with the amputee's impact level.
2. Within the selected column, locate the amputee's weight.

Important Note: If the amputee has a long BK, carries heavy loads or will reach a higher impact level within a year, choose the next category higher.

Warning: Choosing a lower strength category than what is suggested based on the above procedure and patient data will void the warranty. If your patient's weight exceeds the limits of the chart, please browse the catalog for other foot options or call Trulife Customer Service.

	Low	Medium-Low	Medium	Medium-High	High
Category	Walking with Aid	Limited Walking	Walking Uneven Surfaces	Light Sports	Running, Basketball
H9	- -	116-136 kg 256-300 lb	114-136 kg 251-300 lb	91-106 kg 201-235 lb	77-90 kg 170-200 lb
H8	- -	87-115 kg 191-255 lb	91-113 kg 201-250 lb	73-90 kg 161-200 lb	62-76 kg 136-169 lb
H7	123-136 kg 271-300 lb	75-86 kg 166-190 lb	69-90 kg 151-200 lb	55-72 kg 121-160 lb	46-61 kg 101-135 lb
H6	82-122 kg 181-270 lb	59-74 kg 131-165 lb	45-68 kg 101-150 lb	36-54 kg 81-120 lb	29-45 kg 66-100 lb
H5	< 82 kg < 181 lb	< 59 kg < 131 lb	< 45 kg < 101 lb	< 36 kg < 81 lb	< 29 kg < 66 lb



AAA100
Adjustable Height
Clamp Adapter
(page 16)

AAA135
Variable Height
Adapter
(page 17)



Seattle Kinetic Edge with Pylon

- All the benefits of the Kinetic Light with the addition of a removable, sandal toe foot shell and a more lively keel
- Low profile composite keel with multi-axial ankle
- Male pyramid attachment
- Ideal for patients on the edge of advancing to a K3 product

Indications: K2 impact level, maximum weight of 160 kg / 350 lb*

Selection

Model No.	Description	Size	Weight Limit
SKF300	Seattle Kinetic Edge, Light Foot Shell	22-30 cm	160 kg / 350 lb*
SKF303	Seattle Kinetic Edge, Dark Foot Shell	22-30 cm	160 kg / 350 lb*
SKC300	Replacement Foot Shell, Light	22-30 cm	N/A
SKC303	Replacement Foot Shell, Dark	22-30 cm	N/A

Specifications

Weight Limit	160 kg / 350 lb*
Product Sizes	22-30 cm
Toe	Sandal Toe
Foot Shell	Light, Dark
Product Weight	590 g / 1.3 lb (Size 27 without Pylon)
Min. Clearance	7 cm / 2.9 in
Heel Rise	9.5 mm / 3/8 in
Warranty Foot: Foot Shell:	2 years 6 months
Includes	Female Pylon (AAA200)
Made in USA	

*For category #5

Bumper Selection Instructions

To optimize selection, follow the three steps below to determine the appropriate category.

1. Locate the column associated with the amputee's foot length.
2. Locate the row corresponding to the amputee's weight.
3. The area where the column and row intersect lists the appropriate category.

Important Note: If the amputee has a long BK or carries heavy loads choose the next category higher.

Warning: Choosing a lower category than what is suggested based on the above procedure and patient data will void the warranty. If your patient's weight exceeds the limits of the chart, please browse the catalog for other foot options or call Trulife Customer Service.

Foot Length	22	23	24	25	26	27	28	29	30
Weight									
137-160 kg 301-350 lb	4	4	5	5	5	5	5	5	5
101-136 kg 221-300 lb	3	3	4	4	4	4	4	4	4
81-100 kg 177-220 lb	2	2	3	3	3	3	3	3	3
61-80 kg 133-176 lb	1	1	2	2	2	2	2	2	2
< 60 kg < 132 lb	1	1	1	1	1	1	1	1	1



AAA225
Clamp Adapter
(page 16)





AAA100
Adjustable Height
Clamp Adapter
(page 16)

Seattle Kinetic Light with Pylon

- K2 foot with integrated multi-axial ankle
- Proximal male pyramid
- Full length, composite keel
- Sandal toe foot shell (non-removable)
- Integrated ankle function for all categories is the same as the existing Seattle Kinetic

Indications: Moderate impact level, K2, maximum weight of 160 kg /350 lb*

Selection

Model No.	Description	Size	Weight Limit
SKF100	 Seattle Kinetic Light, Light Foot Shell	22-30 cm	160 kg / 350 lb*
SKF103	 Seattle Kinetic Light, Dark Foot Shell	22-30 cm	160 kg / 350 lb*

Bumper Selection Instructions

To optimize selection, follow the three steps below to determine the appropriate category.

1. Locate the column associated with the amputee's foot length.
2. Locate the row corresponding to the amputee's weight.
3. The area where the column and row intersect lists the appropriate category.

Important Note: If the amputee has a long BK or carries heavy loads choose the next category higher.

Warning: Choosing a lower category than what is suggested based on the above procedure and patient data will void the warranty. If your patient's weight exceeds the limits of the chart, please contact Trulife Customer Service.

Foot Length	22	23	24	25	26	27	28	29	30
Weight									
137-160 kg 301-350 lb	4	4	5	5	5	5	5	5	5
101-136 kg 221-300 lb	3	3	4	4	4	4	4	4	4
81-100 kg 177-220 lb	2	2	3	3	3	3	3	3	3
61-80 kg 133-176 lb	1	1	2	2	2	2	2	2	2
< 60 kg < 132 lb	1	1	1	1	1	1	1	1	1



Specifications

Weight Limit	160 kg / 350 lb*
Product Sizes	22-30 cm
Toe	Sandal Toe
Foot Shell	Light, Dark
Product Weight Cat 1-2: Cat 3-5:	481 g / 1.1 lb 564 g / 1.2 lb (Size 25 without Pylon)
Min. Clearance	7.3 cm / 2.9 in
Heel Rise	9.5 mm / 3.8 in
Warranty Foot: Foot Shell:	2 years 1 year
Includes	Female Pylon (AAA200)
Made in USA	

*For category #5



AAA225
Clamp Adapter
(page 16)

AAA100
Adjustable Height
Clamp Adapter
(page 16)





Seattle Kinetic with Pylon

- K2 impact level
- Integrated multi-axial ankle
- Composite keel
- Sandal toe foot shell (non-removable)
- 43 cm / 17 inch male bonded pylon and female foot adapter included
- 160 kg / 350 lb* weight limit

Indications: Moderate impact level, K2; ideal for single speed, community ambulators

Selection

Model No.	Description	Size	Weight Limit
SKF200	 Seattle Kinetic, Light Foot Shell	22-30 cm	160 kg / 350 lb*
SKF203	 Seattle Kinetic, Dark Foot Shell	22-30 cm	160 kg / 350 lb*

Bumper Selection Instructions

To optimize selection, follow the three steps below to determine the appropriate category.

1. Locate the column associated with the amputee's foot length.
2. Locate the row corresponding to the amputee's weight.
3. The area where the column and row intersect lists the appropriate category.

Important Note: If the amputee has a long BK or carries heavy loads choose the next category higher.

Warning: Choosing a lower category than what is suggested based on the above procedure and patient data will void the warranty. If your patient's weight exceeds the limits of the chart, please contact Trulife Customer Service.

Foot Length	22	23	24	25	26	27	28	29	30
Weight									
137-160 kg 301-350 lb	4	4	5	5	5	5	5	5	5
101-136 kg 221-300 lb	3	3	4	4	4	4	4	4	4
81-100 kg 177-220 lb	2	2	3	3	3	3	3	3	3
61-80 kg 133-176 lb	1	1	2	2	2	2	2	2	2
< 60 kg < 132 lb	1	1	1	1	1	1	1	1	1

Specifications

Weight Limit	160 kg / 350 lb*
Product Sizes	22-30 cm
Toe	Sandal Toe
Foot Shell	Light, Dark
Product Weight	495 g / 1.1 lb (Size 25 without Pylon)
Min. Clearance	7 cm / 2.8 in
Heel Rise	9.5 mm / 3/8 in
Warranty Keel:	2 years
Foot Shell:	1 year
Includes	Male Bonded Pylon AAA202
Made in USA	

*For category #5



AAA225
Clamp Adapter
(page 16)

AAA100
Adjustable Height
Clamp Adapter
(page 16)



Seattle Natural with Pyramid

- Lightweight
- One standard Delrin keel per foot size
- Propels amputee through cadence
- Natural looking Lightfoot Foot Shell and realistically detailed Lifecast Foot Shell available in both light and dark (non-removable)

Indications: Low to medium-low impact levels, K1-K2, maximum weight of 102 kg / 225 lb; ideal for single speed, home ambulators or new amputees

Specifications

Weight Limit	102 kg / 225 lb
Product Sizes	
Lightfoot Foot Shell:	22-30 cm
Lifecast Foot Shell:	22-29 cm
Toe	Sandal Toe
Foot Shell	Light, Dark
Product Weight	448 g / 1 lb (Size 25)
Min. Clearance	8.5 cm / 3.3 in
Heel Rise	
SNF160, 163:	9.5 mm / 3/8 in
SNF170, 173:	19 mm / 3/4 in
Warranty	
Keel:	2 years
Foot Shell:	1 year
Pre-assembled with	M10 Foot Bolt and Foot Pyramid
Made in USA	

Selection

Model No.	Description	Size	Weight Limit
SNF160	Seattle Natural with Pyramid / 3/8 in Heel, Lightfoot Foot Shell, Light	22-30 cm	102 kg / 225 lb
SNF163	Seattle Natural with Pyramid / 3/8 in Heel, Lightfoot Foot Shell, Dark	22-30 cm	102 kg / 225 lb
SNF170	Seattle Natural with Pyramid / 3/4 in Heel, Lifecast Foot Shell, Light	22-29 cm	102 kg / 225 lb
SNF173	Seattle Natural with Pyramid / 3/4 in Heel, Lifecast Foot Shell, Dark	22-29 cm	102 kg / 225 lb



Lightfoot Foot Shell,
Light



Lightfoot Foot Shell,
Dark



Lifecast Foot Shell,
Light



Lifecast Foot Shell,
Dark



AAA100
Adjustable Height
Clamp Adapter
(page 16)

AAA200
Aluminum Pylon
Bonded to Female
Adapter (page 17)



Seattle SACH with Pyramid

- Easy rollover
- Lightweight
- Natural looking Lightfoot Foot Shell and realistically detailed Lifecast Foot Shell available in both light and dark (non-removable)

Indications: Low to medium-low impact level, K1, maximum weight 102 kg / 225 lb; ideal for new amputees

Selection

Model No.	Description	Size	Weight Limit
SCH110	 Seattle SACH with Pyramid / 3/8 in Heel, Lightfoot Foot Shell, Light	22-30 cm	102 kg / 225 lb
SCH113	 Seattle SACH with Pyramid / 3/8 in Heel, Lightfoot Foot Shell, Dark	22-30 cm	102 kg / 225 lb
SCH120	 Seattle SACH with Pyramid / 3/4 in Heel, Lifecast Foot Shell, Light	22-29 cm	102 kg / 225 lb
SCH123	 Seattle SACH with Pyramid / 3/4 in Heel, Lifecast Foot Shell, Dark	22-29 cm	102 kg / 225 lb

Specifications

Weight Limit	102 kg / 225 lb
Product Sizes	
Lightfoot Foot Shell:	22-30 cm
Lifecast Foot Shell:	22-29 cm
Toe	Sandal Toe
Foot Shell	Light, Dark
Product Weight	549 g / 1.2 lb (Size 25)
Min. Clearance	
SCH110, 113:	8.5 cm / 3.3 in
SCH120, 123 (Size 22-26):	9 cm / 3.5 in
SCH120, 123 (Size 27-29):	9.5 cm / 3.7 in
Heel Rise	
SCH110, 113:	9.5 mm / 3/8 in
SCH120, 123:	19 mm / 3/4 in
Warranty	
Keel:	2 years
Foot Shell:	1 year
Pre-assembled with	M10 Foot Bolt and Foot Pyramid
Made in USA	



Lightfoot Foot Shell, Light



Lightfoot Foot Shell, Dark



Lifecast Foot Shell, Light



Lifecast Foot Shell, Dark



AAA100
Adjustable Height
Clamp Adapter
(page 16)

AAA200
Aluminum Pylon
Bonded to Female
Adapter (page 17)



Seattle Select Stance Flexion Knee

- Polycentric, pneumatic 5-bar knee
- 9° of adjustable stance flexion
- Shock absorption at heel strike
- Independently adjustable pneumatic flexion and extension control
- Available with built-in proximal pyramid or threaded adapter
- Equipped with a built-in knee cap, preventing clothing from getting caught in the mechanism

Indications:

SSK615: Above knee amputation, medium to high impact level, K3-K4, maximum weight of 100 kg / 220 lb

SSK615-THR: Above knee amputation, medium to high impact level, K3-K4, maximum weight of 125 kg / 275 lb

Selection

Model No.	Description	Weight Limit
SSK615	Seattle Select Stance Flexion Knee with Proximal Pyramid	100 kg / 220 lb
SSK615-THR	Seattle Select Stance Flexion Knee with Threaded Proximal Adapter	125 kg / 275 lb
SSK615-CAP	Knee Cap for Seattle Select Stance Flexion Knees	

Specifications

SSK615	
Weight Limit	100 kg / 220 lb
Product Weight	1100 g / 2.5 lb
Overall Height	24 cm / 9.4 in
Distance to Knee Center	3 cm / 1.2 in
Flexion	145°
Warranty	2 years
SSK615-THR	
Weight Limit	125 kg / 275 lb
Product Weight	1100 g / 2.5 lb
Overall Height	23.2 cm / 9.1 in
Distance to Knee Center	2 cm / 0.8 in
Flexion	145°
Warranty	2 years

Trulife Complementary Items



AAASS237
3-Prong Rotatable
Lamination Piece
(page 16)

AAA214-01
Rotatable Female
Insert
(page 16)



STF190, 193
Seattle Triumph
(page 4)



AAA135
Variable Height
Adapter
(page 17)



Seattle Natural Stance Flexion Knee

- 7° of adjustable stance flexion
- Friction swing phase control
- Shock absorption at heel strike
- Extension assist
- Built-in proximal pyramid and distal tube adapter

Indications: Above knee amputation, low to moderate impact level, K1-K2, maximum weight of 100 kg / 220 lb

Selection

Model No.	Description	Weight Limit
SSK617	Seattle Natural Stance Flexion Knee	100 kg / 220 lb

Specifications

Weight Limit	100 kg / 220 lb
Product Weight	783 g / 1.7 lb
Overall Height	18.5 cm / 7.3 in
Distance to Knee Center	3.2 cm / 1.3 in
Flexion	150°
Warranty	2 years



AAASS237
3-Prong Rotatable
Lamination Piece
(page 16)

AAA214-01
Rotatable Female
Insert
(page 16)

SNF160, 163, 170, 173
Seattle Natural
with Pyramid
(page 12)

AAA200
Aluminum Pylon
Bonded to Female
Adapter (page 17)



Seattle PC Friction Control Knee with Lock

- Polycentric, 4-bar friction control knee with manual lock
- A lightweight, durable knee that can be used in locked position or free motion
- Constant friction braking mechanisms and spring-loaded extension assist helps to create a natural gait
- Built-in proximal pyramid

Indications: Above knee amputation, low to moderate impact level, K1-K2, maximum weight of 136 kg / 300 lb

Selection

Model No.	Description	Weight Limit
SSK602L	PC Friction Control Knee with Lock	136 kg / 300 lb

Specifications

Weight Limit	136 kg / 300 lb
Product Weight	630 g / 1.4 lb
Overall Height	13.5 cm / 5.3 in
Distance to Knee Center	3.2 cm / 1.3 in
Flexion	145°
Warranty	2 years



AAASS237
3-Prong Rotatable
Lamination Piece
(page 16)

AAA214-01
Rotatable Female
Insert
(page 16)

SKF200,300
Seattle Kinetic
with Pylon
(pages 9, 11)

AAA 3-Prong Adapter



Model No.	Description	Rotatable	Material	Product Weight	Weight Limit
AAASS237	3-Prong Laminating Piece SS	Yes	Stainless Steel	107 g / 3.8 oz	136 kg / 300 lb

AAA 4-Hole Adapters



Model No.	Description	Material	Connection	Hole Inner Diameter	Product Weight	Weight Limit
AAA210	4-Hole Pyramid Adapter	Aluminum	Male	10 mm / 0.4 in	90 g / 3.2 oz	160 kg / 350 lb
AAA211	4-Hole Female Adapter**	Aluminum	Female	N/A	88 g / 3.1 oz	160 kg / 350 lb
AAA212	4-Hole Rotatable Flange**	Aluminum	N/A	N/A	58 g / 2.0 oz	160 kg / 350 lb
AAA214-01	Rotatable Female Insert*	Aluminum	Female	N/A	67 g / 2.3 oz	160 kg / 350 lb
AAA216-01	Rotatable Male Pyramid Insert	Aluminum	Male	10 mm / 0.4 in	54 g / 1.9 oz	160 kg / 350 lb
AAA208	4-Hole Double Offset Pyramid Adapter**	Aluminum	Male	10 mm / 0.4 in	90 g / 3.2 oz	160 kg / 350 lb
AAA209	4-Hole Single Offset Pyramid Adapter**	Aluminum	Male	10 mm / 0.4 in	90 g / 3.2 oz	160 kg / 350 lb

*Includes set screws. **Includes mounting screws.

AAA Clamp Adapters



Model No.	Description	Inner Diameter	Material	Connection	Product Weight	Weight Limit
AAA100	Adjustable Height Clamp Adapter	30 mm	Aluminum	Female	120 g / 4.2 oz	160 kg / 350 lb
AAA225	Clamp Adapter	30 mm	Aluminum	Female	90 g / 3.2 oz	160 kg / 350 lb

Made in USA

AAA Double-Ended Adapters



Model No.	Description	Height	Material	Connection	Product Weight	Weight Limit
AAA221	Double-Ended Adapter	35 mm / 1.375 in	Aluminum	Female	109 g / 3.8 oz	160 kg / 350 lb
AAA222	Double-Ended Adapter	45 mm / 1.8 in	Aluminum	Female	125 g / 4.4 oz	160 kg / 350 lb
AAA223	Double-Ended Adapter	60 mm / 2.4 in	Aluminum	Female	142 g / 5 oz	160 kg / 350 lb
AAA224	Double-Ended Adapter	75 mm / 3 in	Aluminum	Female	158 g / 5.6 oz	160 kg / 350 lb

AAA Foot Pyramid



Model No.	Description	Material	Connection	No. of Grooves	A-P Length	Product Weight	Weight Limit
AAA285	Foot Pyramid	Aluminum	Male M10 Thread	One (1)	51 mm / 2 in	67 g / 2.4 oz	136 kg / 300 lb

AAA MODULAR COMPONENTS

AAA Seattle Mounting Plate



Model No.	Description	Material	Outer Diameter	Center Hole	Product Height	Product Weight	Weight Limit
AAASMP (Also on page 22)	Seattle Mounting Plate*	Aluminum	6.9 cm / 2.7 in	M18 Threaded	17 mm / 0.68 in	118 g / 4.2 oz	205 kg / 450 lb

*Includes Expulsion Valve Insert.

AAA Pylons




Model No.	Description	Adapter Material	Connection	Diameter	Length	Product Weight	Weight Limit
AAA200	Aluminum Pylon Bonded to Female Adapter*	Aluminum	Female	30 mm	42 cm / 16.5 in	324 g / 11.4 oz	160 kg / 350 lb
AAA202	Aluminum Pylon Bonded to Male Adapter	Aluminum	Male	30 mm	42 cm / 16.5 in	329 g / 11.6 oz	160 kg / 350 lb
AAA135	Variable Height Adapter*	Aluminum	Female	30 mm	178 mm / 7 in	190 g / 6.7 oz	160 kg / 350 lb

Made in USA

*Includes set screws.

AAAHDTi Heavy-Duty Adapters

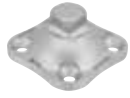

Model No.	Description	Material	Connection	Hole Inner Diameter	Product Weight	Weight Limit
AAAHDTi210	4-Hole Pyramid Adapter	Titanium	Male	N/A	100 g / 3.5 oz	205 kg / 450 lb
AAAHDTi217	Splined 4-Hole Rotatable Adapter	Titanium	Male	N/A	128 g / 4.5 oz	205 kg / 450 lb
AAAHDTi218	Splined 4-Hole Rotatable Adapter	Titanium	Female	N/A	138 g / 4.9 oz	205 kg / 450 lb



Model No.	Description	Inner Diameter	Material	Connection	Product Weight	Weight Limit
AAAHDTi225	Adjustable Height Clamp Adapter	30 mm	Titanium	Female	102 g / 3.6 oz	204 kg / 450 lb



Model No.	Description	Height	Material	Connection	Product Weight	Weight Limit
AAAHDTi221	Double-Ended Adapter	35 mm / 1.38 in	Titanium	Female	118 g / 4.2 oz	204 kg / 450 lb



Model No.	Description	Adapter Material	Connection	Diameter	Length	Product Weight	Weight Limit
AAAHDTi200	Bonded Pylon 16 Inch	Titanium/Composite	Female	30 mm	42 cm / 16.5 in	328 g / 11.6 oz	205 kg / 450 lb





Green Dot Total Contact Suction Socket Valve

- Reliable seal
- Lightweight design

Indications: Transfemoral amputees who use thermoformed sockets

Selection

Model No.	Description
P12-310-1000	Green Dot Total Contact Suction Socket Valve, Plastic, Manual Expulsion



Green Dot SFS Valve

- Lightweight, replaceable valve housing
- Reliable seal
- Lightweight design

Indications: Transfemoral amputees who use thermoformed sockets

Selection

Model No.	Description
P12-340-1000	Green Dot SFS Valve, Plastic, Manual Expulsion



Green Dot Total Contact Suction Valve

- Industry standard in reliable suction
- Lightweight design

Indications: Transfemoral amputees who use lamination sockets

Selection

Model No.	Description
P12-320-1000	Green Dot Total Contact Suction Valve, Stainless Steel, Manual Expulsion

Specifications

(All Green Dot Valves)

Distal End of Valve to Top of Threads	10.8 mm / 0.4 in
Diameter at Valve Base	22 mm / 0.9 in
Warranty	2 years
Includes	Valve Key, Valve Plug, Flat Washer, Screw
Made in USA	



Below Knee Valve

- Self expulsion valve
- Installed at the distal end of socket
- For use with suction suspension

Indications: Transtibial amputees

Selection

Model No.	Description
P36110	Below Knee Valve

Specifications

Minimum Distal Socket Thickness	19 mm / 0.8 in
Warranty	2 years
Made in USA	

OTHER



Specifications

Product Weight	
SEC201, 210:	246 g / 8.6 oz (size 26)
SEC202, 211:	348 g / 12.2 oz (size 26)
SEC220, 223:	94 g / 3.3 oz (size 15)
Product Height	
SEC201, 210:	26 cm / 10.3 in
SEC202, 211:	39 cm / 15.5 in
SEC220, 223:	12-30 cm / 4.8-12 in (varies with foot size)
Warranty	2 years
Made in USA	

Seattle Endoskeletal Cosmesis

- Universal cosmesis shape for temporary or permanent application
- Made of durable urethane for quick, easy finishing
- Easy access for alignment changes
- Please have side and foot size available when ordering

Indications: Transtibial amputations; for use with the following Seattle feet: Lightfoot with Pyramid, Natural with Pyramid, SACH with Pyramid or Child's Play Energy

Selection

Model No.	Description	Foot Size
SEC201	Endoskeletal Cosmesis, Light	22-26 cm
SEC210	Endoskeletal Cosmesis, Dark	22-26 cm
SEC202	Endoskeletal Cosmesis, Light	26-29 cm
SEC211	Endoskeletal Cosmesis, Dark	26-29 cm
SEC220	Child's Play Endoskeletal Cosmesis, Light	13-21 cm
SEC223	Child's Play Endoskeletal Cosmesis, Dark	13-21 cm



Specifications

Product Weight	
SAB310:	426 g / 15.0 oz
SAB320:	369 g / 13.0 oz
Product Height	
SAB310:	12.4 cm / 4.9 in
SAB320:	8.6 cm / 3.4 in
Max. Clearance (including Foot)	15.5 cm / 6.1 in
Core Diameter	5.1 cm / 2.0 in
Largest Width	11.4 cm / 4.5 in
Warranty	2 years
Made in USA	

Seattle Ankle Block

- Lightweight and durable to withstand the forces imposed by dynamic response feet
- High strength, rigid polyurethane foam surrounding a nylon core

Indications:

SAB310: For use with the Seattle LightFoot² models SLF195 and SL198

SAB320: For use with the Seattle Natural Foot models SNF160, SNF163, SNF170, SNF173 and Seattle SACH Foot models SCH110, SCH113, SCH120, SCH123

Selection

Model No.	Description	Weight Limit
SAB310	Seattle Ankle Block, Standard	136 kg / 300 lb
SAB320	Seattle Ankle Block, Low Profile	102 kg / 225 lb



Specifications

Product Weight	
SSY300:	42 g / 1.4 oz
SSY301:	14 g / 0.4 oz
SSY302:	20 g / 0.7 oz
SSY303:	28 g / 0.9 oz
Product Height	
SSY300:	18 mm / 0.7 in
SSY 301, 302, 303:	13 mm / 0.5 in
Min. Clearance	(Including Foot)
SSY300:	5 cm / 2 in
SSY301, 302, 303:	4.4 cm / 1.7 in
Warranty	2 years
Made in USA	

Seattle Symes Nut

- Provides a broad, solid anterior support surface for laminated sockets
- Fabricated from aluminum alloy
- Abrasive etch finish ensures secure socket adhesion during fabrication

Indications:

SSY300: for use with Seattle Lightfoot with Pyramid

SSY301: for use with Child's Play Energy, size 13-15 cm

SSY302: for use with Child's Play Energy, size 16-18 cm

SSY303: for use with Child's Play Energy, size 19-21 cm

Selection

Model No.	Description	Weight Limit
SSY300	Seattle Symes Nut	136 kg / 300 lb
SSY301	Child's Play Symes Nut, Foot Size 13-15 cm	31 kg / 69 lb
SSY302	Child's Play Symes Nut, Foot Size 16-18 cm	40 kg / 89 lb
SSY303	Child's Play Symes Nut, Foot Size 19-21 cm	65 kg / 144 lb



Specifications

Product Weight	118 g / 4.2 oz
Product Height	17 mm / 0.68 in
Outer Diameter	6.9 cm / 2.7 in
Center Hole	M18 Threaded
Mounting Pattern	4 x M6 Threaded
Bolt Circle	50.8 mm / 2.0 in
Warranty	2 years
Made in USA	

AAA Seattle Mounting Plate

- Fabricated from aluminum alloy
- Includes expulsion valve assembly
- Lamination plug (included) keeps the threaded center hole free of resin during the lamination process, and serves as a guide for drilling the final socket vent hole
- Lamination dam keeps the distal surface clean and free of resin during the lamination process
- Expulsion valve and lamination plug do not require any special tool

Selection

Model No.	Description	Weight Limit
AAASMP	Seattle Mounting Plate	205 kg / 450 lb



Note: Foot shell is not removable.

Specifications

Weight Limit	65 kg / 144 lb*
Product Sizes	13-21 cm
Toe	13-15 cm: Solid Toe 16-21 cm: Sandal Toe
Foot Shell	Light, Dark
Product Weight	150 g / 5 oz (Size 17)
Min. Clearance	
Size 13-15:	4 cm / 1.5 in
Size 16-18:	4.5 cm / 1.8 in
Size 19-21:	4.5 cm / 1.8 in
Heel Rise	
Size 13-15:	6.3 mm / 1/4 in
Size 16-18:	9.5 mm / 3/8 in
Size 19-21:	12.7 mm / 1/2 in
Warranty	
Keel:	2 years
Foot Shell:	1 year
Includes	M8 x 80 mm Foot Bolt Kit (SFB546)
Made in USA	

*For #6 Keel

Seattle Child's Play Energy

- Pediatric model of the Seattle Energy
- Dynamic toe and cushioned heel response
- M8 foot bolt attachment
- Non-removable, sandal-toe foot shell (size 16-21 cm)
- High weight limit of 65 kg / 144 lb

Indications: High impact level, K4, maximum weight of 65 kg / 144 lb*

Selection

Model No.	Description	Size	Weight Limit
SEF132	Child's Play Energy, Light Foot Shell	13-21 cm	65 kg / 144 lb*
SEF144	Child's Play Energy, Dark Foot Shell	13-21 cm	65 kg / 144 lb*

Keel Selection Instructions

To determine which keel is appropriate for your patient, check the below table for keel availability by foot size. Please note that sizes 13-15 feet do not have sandal toe foot shells.

Important Note: Choose the next higher keel level if your patient has a long BK limb, regularly carries heavy loads, or is at the upper limit for weight or impact level.

Warning: Choosing a lower strength category than what is suggested based on the above procedure and patient data will void the warranty. If your patient's weight exceeds the limits of the chart please call Trulife Customer Service or your Territory Manager.

Foot Length	13	14	15	16	17	18	19	20	21
Weight Limit									
Keel									
< 15 kg < 33 lb	1	✓	✓	✓					
16-25 kg 34-55 lb	2	✓	✓	✓	✓	✓	✓		
26-35 kg 56-77 lb	3	✓	✓	✓	✓	✓	✓	✓	✓
36-45 kg 78-99 lb	4				✓	✓	✓	✓	✓
46-55 kg 100-121 lb	5						✓	✓	✓
56-65 kg 122-144 lb	6						✓	✓	✓



AAAK285
Child's Play Energy
Foot Pyramid
(page 24)

AAAK100
Child's Play
Adjustable Height
Clamp Adapter
(page 24)

AAAK135
Child's Play Variable
Height Adapter
(page 24)

SSK610A
Child's Play Knee
with Manual Lock
(page 23)



Seattle Child's Play Knee with Manual Lock

- Lightweight, 4-bar, polycentric design
- Friction swing control and extension assist
- Adjustable stability control
- Free motion setting available
- Rotatable knee disarticulation adapter available (SSK610A-DIS)

Indications: AK amputation, pediatric amputee, maximum weight of 55 kg / 121 lb

Selection

Model No.	Description	Weight Limit
SSK610A	Child's Play Knee with Manual Lock	55 kg / 121 lb
SSK610A-DIS	Child's Play Lamination Bracket	

Specifications

Weight Limit	55 kg / 121 lb
Product Weight	310 g / 11 oz
Overall Height	12 cm / 4.7 in
Distance to Knee Center	3 cm / 1.2 in
Flexion	140°
Warranty	2 years
Includes	Pyramid Adapter

Replacement Parts

Model No.	Description
SSK610A-KIT	Refurbishment Kit
SSK610A-PYR	Replacement Pyramid
SSK610A-BMP	Replacement Bumper Kit



AAAK285
Child's Play
Energy
Foot Pyramid
(page 24)







AAAK221
Child's Play
Double-Ended
Adapter
(page 24)

SSK610A-DIS
Child's Play
Lamination
Bracket


AAAK135
Child's Play
Variable
Height Adapter
(page 24)

SEF132, 144
Child's Play
Energy
(page 22)


Child's Play Modular Adapters

	Model No.	Description	Material	Connection	Diameter	Length	Product Weight	Weight Limit
	AAAK210	Child's Play 4-Hole Pyramid Adapter*	Aluminum	Male	N/A	N/A	49 g / 1.7 oz	65 kg / 144 lb
	AAAK211	Child's Play 4-Hole Female Adapter*	Aluminum	Female	N/A	N/A	39 g / 1.4 oz	65 kg / 144 lb
	AAAK221	Child's Play Double-Ended Adapter	Aluminum	Female	32 mm / 1.26 in	N/A	50 g / 1.8 oz	65 kg / 144 lb
	AAAK100	Child's Play Adjustable Height Clamp Adapter	Aluminum	Female	22 mm (ID)	N/A	58 g / 2 oz	65 kg / 144 lb
	AAAK225	Child's Play Clamp Adapter	Aluminum	Female	22 mm (ID)	N/A	40 g / 1.4 oz	65 kg / 144 lb
	AAAK135	Child's Play Variable Height Adapter	Aluminum	Female	22 mm	178 mm / 7.0 in	87 g / 3.1 oz	65 kg / 144 lb

Child's Play Energy Foot Adapter

	Model No.	Description	Material	Connection	No. of Grooves	Length	Product Weight	Weight Limit
	AAAK285	Child's Play Energy Foot Pyramid	Aluminum	Male M8 threaded	One (1)	37 mm / 1.45 in	33 g / 1.2 oz	65 kg / 144 lb

Child's Play Seattle Mounting Plate

	Model No.	Description	Material	Diameter	Center Hole	Height	Product Weight	Weight Limit
	AAAKSMP	Child's Play Seattle Mounting Plate	Aluminum	5.8 cm / 2.3 in (OD)	31.75 mm / 1.25 in	14 mm / 0.55 in	56 g / 2.0 oz	65 kg / 144 lb

All Child's Play Modular Adapters are made in USA and come with a 2 year warranty.
*Includes mounting screws.



Specifications

Product Weight	94 g / 3.3 oz (Size 15)
Calf Circumference	20 cm / 10.2 in
Height	12-30 cm / 4.8-12.1 in (Height depends on foot size)
Fits Foot Size	13-21 cm
Warranty	2 years
Made in USA	

Child's Play Endoskeletal Cosmesis

- Universal cosmesis shape for temporary or permanent application
- Made of durable urethane for quick, easy finishing
- Easy access for alignment changes
- Please have side and foot size available when ordering

Indications: Transtibial amputations; for use with Child's Play Energy

Selection

Model No.	Description	Foot Size
SEC220	Child's Play Endoskeletal Cosmesis, Light	13-21 cm
SEC223	Child's Play Endoskeletal Cosmesis, Dark	13-21 cm



Specifications

Diameter (OD)	9.2 cm / 3.6 in (OD)
Diameter (Top)	7.1 cm / 2.8 in
Height	3.2 cm / 1.3 in
Product Weight	93 g / 3.2 oz
Warranty	2 years
Made in USA	

Child's Play Laminating Block

- High strength, rigid polyurethane foam around a nylon core
- 4-hole connection
- Easily shaped during fabrication

Indications: For use with Child's Play Energy

Selection

Model No.	Description	Weight Limit
SLB240	Child's Play Laminating Block	65 kg / 144 lb



Specifications

Fits Foot Size	
Size 13-15:	SFS240, 241
Size 16-18:	SFS242, 243
Size 19-21:	SFS244, 245
Warranty	2 years
Made in USA	

Child's Play Foot Spacer

- Provides versatility during setup and fabrication
- Designed of lightweight nylon
- For temporary or permanent use

Indications:

SFS240, 241 – Ideal for use with the Child's Play Energy, size 13-15 cm

SFS242, 243 – Ideal for use with the Child's Play Energy, size 16-18 cm

SFS244, 245 – Ideal for use with the Child's Play Energy, size 19-21 cm

Selection

Model No.	Description	Size	Product Wt	Weight Limit
SFS240	Child's Play Foot Spacer	0.6 cm / 0.25 in	6.0 g / 0.2 oz	31 kg / 69 lb
SFS241	Child's Play Foot Spacer	2.5 cm / 1.0 in	20.0 g / 0.7 oz	31 kg / 69 lb
SFS242	Child's Play Foot Spacer	0.6 cm / 0.25 in	9.0 g / 0.3 oz	41 kg / 89 lb
SFS243	Child's Play Foot Spacer	2.5 cm / 1.0 in	28.0 g / 0.9 oz	41 kg / 89 lb
SFS244	Child's Play Foot Spacer	0.6 cm / 0.25 in	11.0 g / 0.3 oz	65 kg / 144 lb
SFS245	Child's Play Foot Spacer	3.0 cm / 1.2 in	45.0 g / 1.5 oz	65 kg / 144 lb



Specifications

Height	8.6 cm / 3.4 in
Max Clearance	(Including Foot)
SAB314:	9.4 cm / 3.7 in
SAB315:	10.4 cm / 4.1 in
SAB316:	11.2 cm / 4.4 in
Core Diameter	
SAB314:	3.0 cm / 1.2 in
SAB315:	3.5 cm / 1.4 in
SAB316:	4.1 cm / 1.6 in
Largest Width	7.9 cm / 3.1 in
Warranty	2 years
Made in USA	

Child's Play Energy Ankle Block

- Lightweight, durable ankle block
- High strength, rigid polyurethane foam surrounding a nylon core

Indications:

SAB314 – For use with the Child's Play Energy, size 13-15 cm

SAB315 – For use with the Child's Play Energy, size 16-18 cm

SAB316 – For use with the Child's Play Energy, size 19-21 cm

Selection

Model No.	Description	Foot Size	Product Weight	Weight Limit
SAB314	Child's Play Energy Ankle Block	13-15 cm	123 g / 4.3 oz	35 kg / 77 lb
SAB315	Child's Play Energy Ankle Block	16-18 cm	146 g / 5.1 oz	45 kg / 99 lb
SAB316	Child's Play Energy Ankle Block	19-21 cm	179 g / 6.3 oz	65 kg / 144 lb



Specifications

Min Clearance	(Including Foot)
SLC301:	5.6 cm / 2.2 in
SLC302:	5.8 cm / 2.3 in
SLC303:	6.0 cm / 2.4 in
Max Clearance	(Including Foot)
SLC301:	9.4 cm / 3.7 in
SLC302:	10.4 cm / 4.1 in
SLC303:	11.2 cm / 4.4 in
Warranty	2 years
Made in USA	

Child's Play Laminating Core

- High strength, low profile alternative to an ankle block
- Requires minimal keel to socket clearance
- Ribbed nylon for secure adhesion in the lamination process

Indications: For laminated sockets, maximum weight of 65 kg / 144 lb; for use with Child's Play Energy

Selection

Model No.	Description	Foot Size	Product Weight	Weight Limit
SLC301	Child's Play Laminating Core	13-15 cm	28 g / 0.9 oz	31 kg / 69 lb
SLC302	Child's Play Laminating Core	16-18 cm	50 g / 1.7 oz	40 kg / 89 lb
SLC303	Child's Play Laminating Core	19-21 cm	79 g / 2.7 oz	65 kg / 144 lb



Specifications

Min Clearance	(Including Foot)
	4.0 cm / 1.6 in
Warranty	2 years
Made in USA	

Child's Play Symes Nut

- Provides a broad, solid anterior support surface
- Fabricated from aluminum alloy
- Abrasive etch finish ensures secure socket adhesion during fabrication

Indications: For use with Child's Play Energy

Selection

Model No.	Description	Foot Size	Product Weight	Weight Limit
SSY301	Child's Play Symes Nut	13-15 cm	14 g / 0.4 oz	31 kg / 69 lb
SSY302	Child's Play Symes Nu	16-18 cm	120 g / 0.7 oz	40 kg / 89 lb
SSY303	Child's Play Symes Nut	19-21 cm	28 g / 0.9 oz	65 kg / 144 lb



Warranty

Warranty Specifications

Feet: Seattle Triumph Seattle LP	Keel: 3 years Cosmesis: 6 mos
Feet: Seattle Energy Seattle Kinetic Edge	Keel: 2 years Cosmesis: 6 mos
Feet: All Other Models	Keel: 2 years Cosmesis: 1 year
Knees: All Models	2 years
Modular Components: All Models	2 years
Valves: All Models	2 years
Endoskeletal and Exoskeletal Components:	2 years
Repairs and Refurbishes:	90 days

Note: Warranty durations begin the date of installation, repair or refurbish.

Trulife warrants each of its prosthetic devices subject to conditions below, to be free from defects in material and workmanship under normal use, service, and proper installation, for the time period specified below:

If failure of a purchased product occurs during this period because of a defect in material or workmanship, upon satisfaction of the conditions set forth below, Trulife will replace, repair, or refund the product at no charge to the original retail purchaser.

Your sole and exclusive remedy, as the original retail purchaser, in the event of a defect is found, is expressly limited to replacement or repair of the product, although Trulife may elect to refund, not to exceed the purchase price, rather than replace or repair the product. Replacement products are warranted for only the remainder of the original warranty period. Warranty is nontransferable and may only be exercised by the original purchaser.

The limited warranty contained herein is in lieu of all other written or expressed warranties. Any implied warranty, including the implied warranties of merchantability or fitness for a particular purpose, is limited in length to the duration of this limited warranty. Some states do not allow limitation on how long an implied warranty lasts, so the above limitation may not apply to you.

Trulife's sole responsibility is to replace, repair, or refund as stated herein, and Trulife shall not be liable for any special, indirect, or consequential damages including loss of time, inconvenience, loss of the use of the product, or commercial losses. Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation or exclusion may not apply to you.

This warranty gives you specific legal rights, and you may also have other rights which vary from state to state.

The application guidelines for this Trulife product are for the use of and by certified, qualified practitioner only. Patients are not to attempt to apply or adjust the item unless expressly instructed to do so by the practitioner responsible for the prescription and/or initial fitting of the device. All patient questions should be referred to the practitioner and not to the manufacturer. The manufacturer warrants only that the enclosed product has been inspected for quality and can be effective for certain indications, but final decisions and ongoing monitoring must be made by the prosthetic professional(s) prescribing and/or fitting the device to determine its effectiveness for an individual patient. Patient compliance is an integral part of the entire protocol and must be adhered to in order to avoid potential problems and to maximize the effectiveness of the prescribed product.

As a condition of the sale of any Trulife product, this prosthesis is restricted to a "Single Patient Use Only" by the originally fitted patient in order to protect the care provider and the patient against potentially adverse consequences of infectious disease transmission, material instability in adapting to the configuration of the original user and/or decrease in effectivity. Any express or implied warranties are voided if the product is reused or fitted to another patient. Additionally, a license of right to use under any relevant patents pertaining to the product is terminated with the cessation of use by the original patient. As with all Trulife products, this product must be prescribed and applied by a qualified practitioner to determine it meets the needs of the particular patient and accomplishes the desired results.





Trulife United States

2010 East High Street
Jackson, MI 49203
t: 800.492.1088
f: 800.245.3765
info-usa@trulife.com

26296 Twelve Trees Lane NW
Poulsbo WA 98370
t: 360.697.5656
f: 360.697.6843
info-usa@trulife.com

445 Sequoia Drive
Bellingham WA 98226
t: 866.477.6262
f: 360.714.9007
info-usa@trulife.com

Trulife Canada

39 East Davis Street
Trenton, ON K8V 4K8
t: 800.267.2812
f: 613.392.4139
infocanada@trulife.com

Trulife United Kingdom

41 Amos Road, Meadowhall
Sheffield, UK S9 1BX
t: +44 (0) 1142 618100
f: +44 (0) 1142 618188
info@trulife.co.uk

Trulife Ireland

Airton Road
Tallaght
Dublin D24 R859
t: +353 (0) 1 4511755
f: +353 (0) 1 4275790
info@trulife.com

www.trulife.com